

SINGLE USE KIT

STERILE R

NEWCLIP-TECHNICS



INITIAL MTP



Ready
when you are!

With a non sterile
standard kit



Calling on medical staff

Constraints



Complex
traceability



Contracted out
sterilization



Suppliers'
deadline


High costs



\$ Stocks
\$ Control
\$ Cleaning
\$ Decontamination
\$ Sterilization



Bulky storage

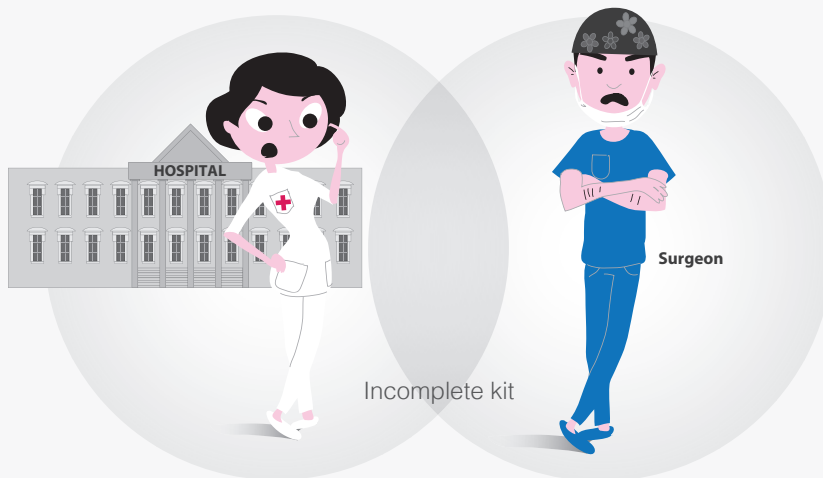

Complex process



Prevents an
effective solution
& a quick response



Defective
sterilization



Incomplete
kit



Damaged
instrumentation



INCREASED
RISKS

NON OPTIMAL
surgery



URGENT SURGICAL CASES
COMPROMISED

With the kit



Safety



Risk
of contamination

Cost efficiency



Controlled stocks
Simplified control


- 0 Cleaning
- 0 Decontamination
- 0 Sterilization



Sundry
expenses



Optimized storage

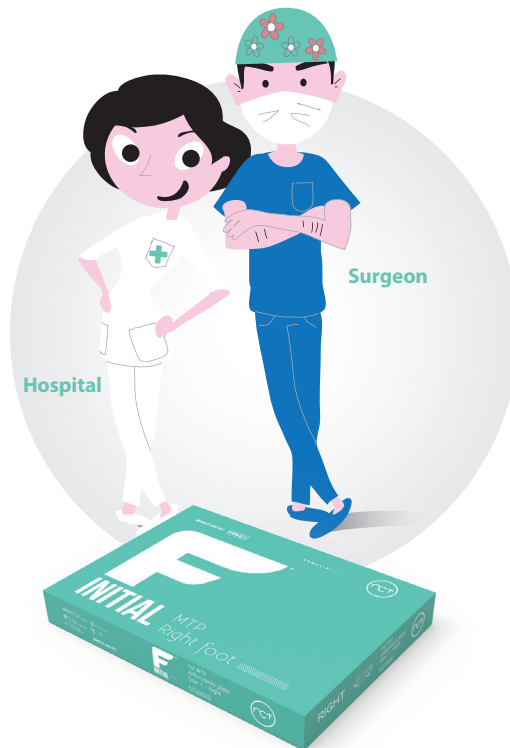

Efficiency



- 1 Delivery
- 2 Storage
- 3 Surgery



An effective
solution &
a quick response



Available when needed



READY-TO-USE FOR
SURGERY



Optimized handling of
URGENT SURGICAL CASES

STERILE R SINGLE USE KIT

with state-of-the-art implants

Ready
when you are!



Safety:

The Initial F - MTP™ kit is fully traceable and has a shelf life of 5 years.
Its instrumentation and implants are “always new” and have never been opened or used before.



innovation



Costs:

Initial F - MTP™ is a cost-effective solution.
The additional costs including cleaning, decontamination, sterilization of kits are cancelled.



Storage:

Initial F - MTP™ kit can be easily stored in the operating room because of its small size.



Buying procedure:

Initial F - MTP™ facilitates buying procedures: restocking and orders are simplified, stock management is optimized.



Contamination:

The combination of sterile implants and sterile single-use instrumentation minimizes contamination risks.

Initial F - MTP™ kits

Technical features

> Indications

The implants of the Initial F - MTP™ range are intended for arthrodeses, fractures and osteotomies fixation and revision surgeries of the foot in adults.

> Contraindications

- Serious vascular deterioration, bone devitalization,
- Pregnancy.
- Acute or chronic local or systemic infections.
- Lack of musculo-cutaneous cover, severe vascular deficiency affecting the concerned area.
- Insufficient bone quality preventing a good fixation of the implants into the bone,
- Muscular deficit, neurological deficiency or behavioral disorders, which could submit the implant to abnormal mechanical strains.
- Allergy to one of the materials used or sensitivity to foreign bodies.
- Serious problems of non-compliance, mental or neurological disorders, failure to follow post-operative care recommendations.
- Unstable physical and/or mental condition.

> Plate for the first metatarso-phalangeal (MTP) joint arthrodesis

Examples of applications: hallux rigidus, severe hallux valgus, polyarthritis

- **Range of precontoured plates:** the design of this implant is the result of a proprietary state-of-the-art mapping technology to establish an optimized congruence between the plate and the bone.




- **Low profile plate:** 1,7 mm thick, thus limiting soft tissue irritation risks while providing an optimized mechanical stability.



- **Hole for pin** to temporarily stabilize the plate.

- **Oblong hole** for pin to achieve compression without removing the pin and to ensure the guiding on the metatarsal.

- **Ramp oblong hole**

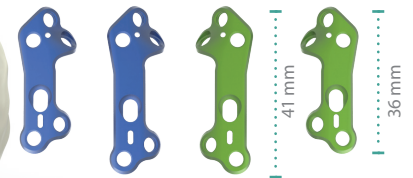

- **2 sizes of plates** for the right (green plates) and left (blue plates) sides offering versatile solutions.



PLATE BENDING

Bending is only possible in the areas intended for this purpose. A bendable area must be bent only once, in one direction and not be performed excessively. The holes must be protected to avoid damaging the fixation.

Initial F - MTP™ kits


Implants - Technical features

> Fixations and screws


- **A single screw diameter:** Ø2.8 mm.
Both locking (SLT2.8Lxx) and non locking screws (RLT2.8Lxx) are available.

- **Screw head is buried in the plate** (1) to limit the risk of soft tissue irritation.

- **The hexalobular screw stamp** improves torque transmission.

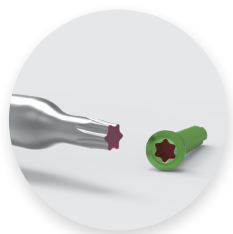


Hexalobular stamp

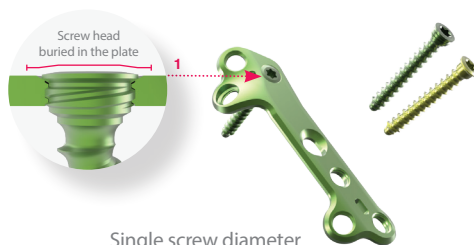

Single screw diameter

> Efficient locking

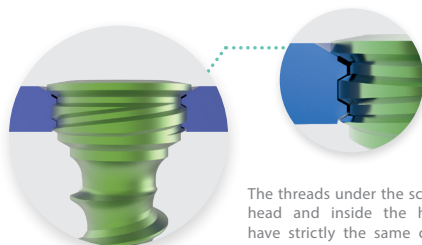

The threads under the screw head and inside the hole have strictly the same characteristics.

Features :


- The screw head is stopped in the hole by its cap, ensuring the locking,
- The screw head is buried in the plate,
- Plate and screws are all made of titanium alloy.

Coaptation of both profiles during locking, limiting a cold welding risk and improving the removal properties.


> Specific fixations for stable assembly

Ramp oblong hole

The ramp oblong hole enables a simple and controlled compression by its screw-plate interface.

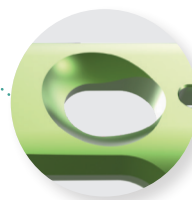


Compression of the joint up to 1.5 mm

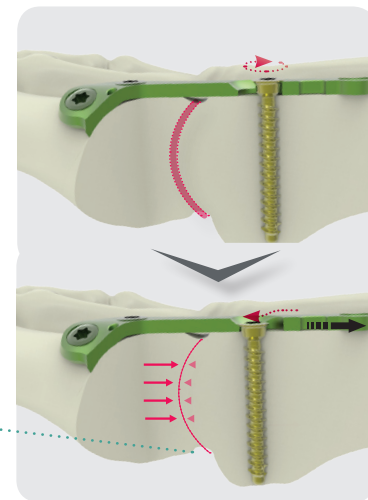


Hole for transfixation screw

The transfixation screw goes through the 1st MTP joint providing stability to the assembly.

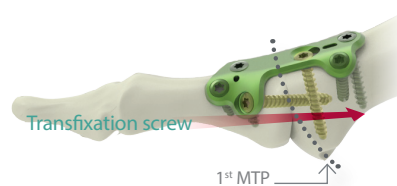

Holes for converging screws in the distal and proximal areas

Converging screws allowing a stable fixation of the system.

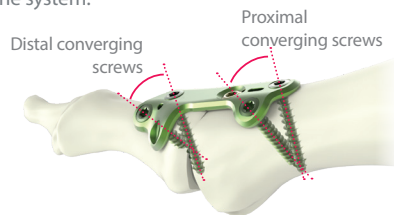


Initial F - MTP™ kits


Instrumentation - Technical features

> Convex and concave reamers

Convex and concave reamers are used respectively to prepare the surfaces of the head of the first metatarsal and the base of the phalanx. Ensuring congruity of the surfaces.

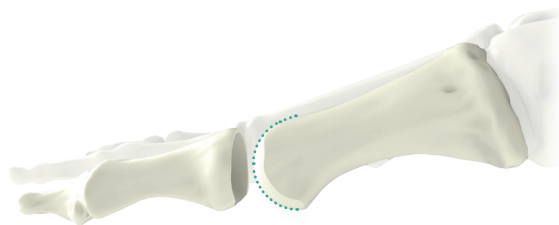


Reamers are available in 3 sizes :

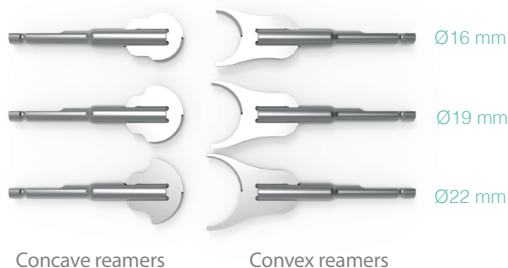

> Handle for guide gauge

Before performing the drilling into the oblong hole, snap the handle for guide gauge on the Ø2.0 mm threaded guide gauge.

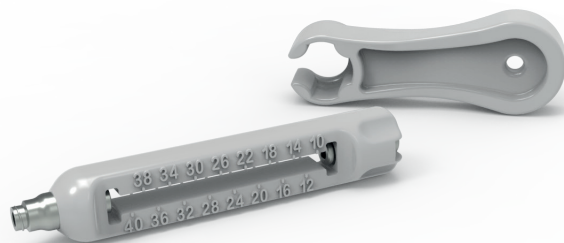

> Templates

The Initial F - MTP™ templates are available separately and allow to quickly and simply determine the appropriate kit.

• REAMER TEMPLATE

The template for **Initial F - MTP™ - Reamers kits** allows to determine the appropriate reamer diameter (Ø16 mm, Ø19mm or Ø22mm) to be used for joint preparation.



• IMPLANT TEMPLATE

The template for **Initial F - MTP kits** allows to determine the desired plate size prior opening a kit.




Initial F - MTP™ kits

Surgical technique - Joint surfaces preparation


Example: surgical technique
with a Ø16 mm reamers kit
(KIT-MI16).



1. Dislocate the joint so as to expose the head of the first metatarsal and the proximal base of the first phalanx.

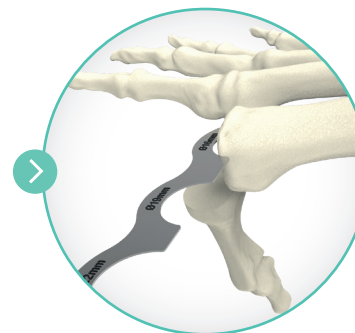


2. Use the reamers template to determine the appropriate reamers kit for joint preparation.

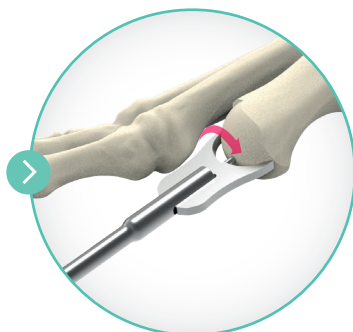

3. Insert the Ø1.6 mm pin through the head of the first metatarsal into the medullary cavity. With the chosen convex reamer, progressively remove the cartilage surface. Then, remove the reamer and the pin.



4. Expose the base of the phalanx and insert the Ø1.6 mm pin so as to achieve the proper alignment with the diaphysis.

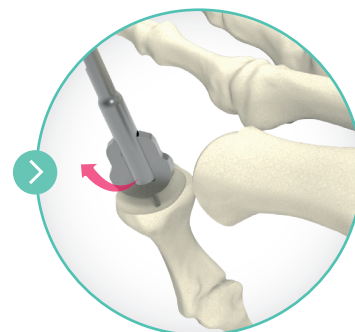

5. Take a concave reamer with **the same diameter** as the convex reamer (determined at step 2). Insert it along the pin and perform the reaming until the cartilage surface has been removed. Then, remove the reamer and the pin.


Initial F - MTP™ kits

Surgical technique


Example: surgical technique with a standard plate, size 1 (KIT-MD1D)




1. Select the correct kit according to the template.

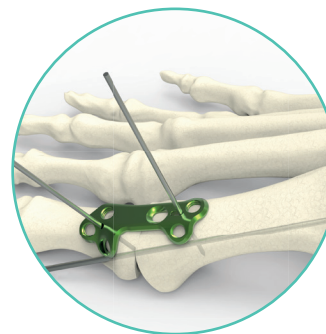


2. Position the joint in the desired direction and stabilize it using a Ø1.6 mm pin. Then, position the plate and stabilize it temporarily by inserting a Ø1.2 mm pin into the dedicated oblong hole pin.

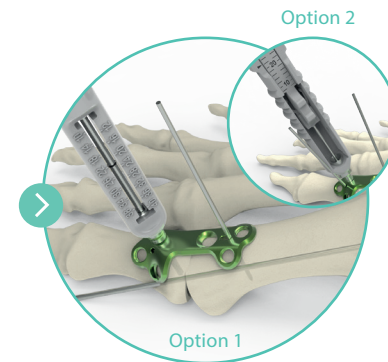


3. Lock the Ø2.0 mm threaded guide gauge into the distal lateral hole and perform the drilling.
Option 1 - Determine the screw length using the gauge.
Option 2 - Determine the screw length using the length gauge.

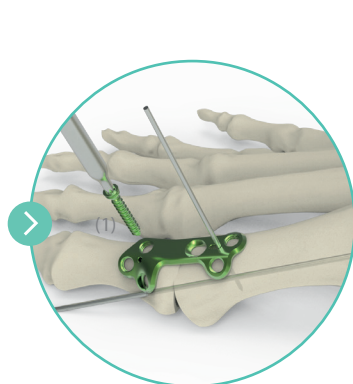

4. Insert the Ø2.8 mm green locking screw with the screwdriver. Repeat the same procedure for the most distal hole (1).

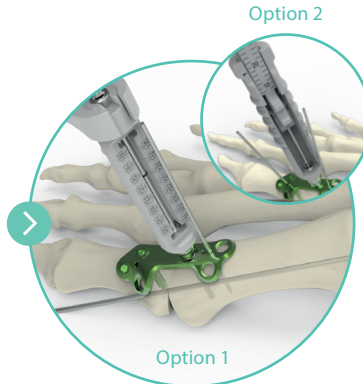

5. Snap the handle for guide gauge and perform the drilling using the assembly into proximal part of the ramp oblong hole.
Option 1 - Determine the screw length using the gauge.
Option 2 - Determine the screw length using the length gauge.

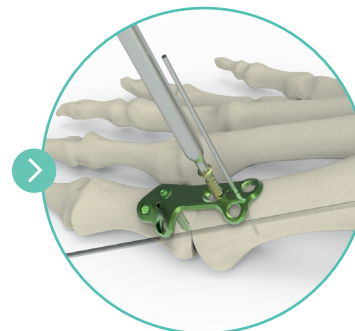

6. In the ramp oblong hole, insert a Ø2.8 mm yellow non locking screw and perform the compression using the screwdriver. Insert the Ø2.8 mm green locking screws into the 2 proximal holes following the steps 2 & 3, then remove the pins.



Final result
 Drill through the distal hole dedicated for the transfixation screw. Finalized the osteosynthesis by inserting a Ø2.8 mm yellow non-locking screw using the screwdriver.

Initial F - MTP™ kits

References



| INITIAL F - MTP™ KITS | |
|-----------------------|--|
| Ref. | Description |
| KIT-MD1D | 1 st MTP Arthrodesis kit - Right - Size 1 |
| KIT-MD1G | 1 st MTP Arthrodesis kit - Left - Size 1 |
| KIT-MD2D | 1 st MTP Arthrodesis kit - Right - Size 2 |
| KIT-MD2G | 1 st MTP Arthrodesis kit - Left - Size 2 |

| INITIAL F - MTP™ KITS - INSTRUMENTATION CONTENT | |
|--|---|
| Description | |
| Ø2.0 mm quick coupling drill bit - L 125 mm | 1 |
| Ø2.0 mm threaded guide gauge for Ø2.8 mm screws | 1 |
| Length gauge for Ø2.8 and Ø3.5 mm screws - L 8-40 mm | 1 |
| T8 prehensor screwdriver | 1 |
| Handle for guide gauge | 1 |
| Pin Ø1.2 L70 mm | 2 |
| Pin Ø1.6 L100 mm | 1 |

NB : Supplemental screws are available in sterile package (cf. : Initial F - MTP™ additional kits, additional implants)

| INITIAL F - MTP™ KITS - IMPLANTS CONTENT | | | QUANTITY PER KIT | | | |
|--|-----------|--|------------------|----------|----------|----------|
| | Ref. | Description | KIT-MD1D | KIT-MD1G | KIT-MD2D | KIT-MD2G |
| PLATES | FMTDD1 | 1 st MTP arthrodesis plate - Size 1 - Right | 1 | - | - | - |
| | FMTGD1 | 1 st MTP arthrodesis plate - Size 1 - Left | - | 1 | - | - |
| | FMTDD2 | 1 st MTP arthrodesis plate - Size 2 - Right | - | - | 1 | - |
| | FMTGD2 | 1 st MTP arthrodesis plate - Size 2 - Left | - | - | - | 1 |
| LOCKING SCREWS Ø2.8 MM | SLT2.8L14 | Locking screw Ø2.8 mm - L 14 mm | 1 | 1 | 1 | 1 |
| | SLT2.8L16 | Locking screw Ø2.8 mm - L 16 mm | 2 | 2 | 2 | 2 |
| | SLT2.8L18 | Locking screw Ø2.8 mm - L 18 mm | 2 | 2 | 2 | 2 |
| | SLT2.8L20 | Locking screw Ø2.8 mm - L 20 mm | 1 | 1 | 1 | 1 |
| NON LOCKING SCREWS Ø2.8 MM | RLT2.8L16 | Non locking screw - Ø2.8 mm - L 16 mm | 1 | 1 | 1 | 1 |
| | RLT2.8L18 | Non locking screw - Ø2.8 mm - L 18 mm | 1 | 1 | 1 | 1 |
| | RLT2.8L22 | Non locking screw - Ø2.8 mm - L 22 mm | 1 | 1 | 1 | 1 |
| | RLT2.8L26 | Non locking screw - Ø2.8 mm - L 26 mm | 1 | 1 | 1 | 1 |

Initial F - MTP™ - Additional sterile implants and kits

References

Additional implants

Sterile screws packaged in the Supplemental sterile screw caddy



| NON LOCKING SCREWS - Ø2.8 mm* | | | LOCKING SCREWS - Ø2.8 mm* | | |
|-------------------------------|---|-----|---------------------------|---|-----|
| Ref. | Description | Qty | Ref. | Description | Qty |
| RLT2.8L10-ST | Non locking screw - Ø2.8 mm - L 10 mm - STERILE | 1 | SLT2.8L10-ST | Locking screw - Ø2.8 mm - L 10 mm - STERILE | 1 |
| RLT2.8L12-ST | Non locking screw - Ø2.8 mm - L 12 mm - STERILE | 1 | SLT2.8L12-ST | Locking screw - Ø2.8 mm - L 12 mm - STERILE | 2 |
| RLT2.8L14-ST | Non locking screw - Ø2.8 mm - L 14 mm - STERILE | 1 | SLT2.8L14-ST | Locking screw - Ø2.8 mm - L 14 mm - STERILE | 3 |
| RLT2.8L16-ST | Non locking screw - Ø2.8 mm - L 16 mm - STERILE | 1 | SLT2.8L16-ST | Locking screw - Ø2.8 mm - L 16 mm - STERILE | 3 |
| RLT2.8L18-ST | Non locking screw - Ø2.8 mm - L 18 mm - STERILE | 1 | SLT2.8L18-ST | Locking screw - Ø2.8 mm - L 18 mm - STERILE | 3 |
| RLT2.8L20-ST | Non locking screw - Ø2.8 mm - L 20 mm - STERILE | 1 | SLT2.8L20-ST | Locking screw - Ø2.8 mm - L 20 mm - STERILE | 2 |
| RLT2.8L22-ST | Non locking screw - Ø2.8 mm - L 22 mm - STERILE | 1 | SLT2.8L22-ST | Locking screw - Ø2.8 mm - L 22 mm - STERILE | 2 |
| RLT2.8L24-ST | Non locking screw - Ø2.8 mm - L 24 mm - STERILE | 1 | SLT2.8L24-ST | Locking screw - Ø2.8 mm - L 24 mm - STERILE | 1 |
| RLT2.8L26-ST | Non locking screw - Ø2.8 mm - L 26 mm - STERILE | 1 | SLT2.8L26-ST | Locking screw - Ø2.8 mm - L 26 mm - STERILE | 1 |
| RLT2.8L28-ST | Non locking screw - Ø2.8 mm - L 28 mm - STERILE | 1 | SLT2.8L28-ST | Locking screw - Ø2.8 mm - L 28 mm - STERILE | 1 |
| RLT2.8L30-ST | Non locking screw - Ø2.8 mm - L30 mm - STERILE | 1 | SLT2.8L30-ST | Locking screw - Ø2.8 mm - L30 mm - STERILE | 1 |
| RLT2.8L32-ST | Non locking screw - Ø2.8 mm - L32 mm - STERILE | 1 | SLT2.8L32-ST | Locking screw - Ø2.8 mm - L32 mm - STERILE | 1 |
| RLT2.8L34-ST | Non locking screw - Ø2.8 mm - L34 mm - STERILE | 1 | SLT2.8L34-ST | Locking screw - Ø2.8 mm - L34 mm - STERILE | 1 |

* Yellow anodized

* Green anodized

Removal and rescue kits


Sterile instruments

| REMOVAL AND RESCUE KITS | | |
|-------------------------|-----------------------------|--|
| Ref. | Description | Content |
| KIT-REMOVE-2 | Removal kit for T8 hexalobe | - 1 x T8 Prehensor Screwdriver |
| KIT-RESCUE-4 | Rescue kit for Initial MTP | - 1x Ø2.0 mm threaded guide gauge for Ø2.8 mm screws - 1x Ø2.0 mm quick coupling drill bit – L125 mm. |

Additional instrumentation kits

Convex & Concave reamers


| SINGLE USE CONVEX AND CONCAVE REAMERS - STERILE PACKAGING | | |
|---|--|---|
| Ref. | Description | Qty |
| KIT-MI16 | Ø16 mm reamers kit for 1 st MTP arthrodesis | - Ø16 mm Initial convex reamer - Ø16 mm Initial concave reamer - Pin Ø1.6 L100 mm x 2 |
| KIT-MI19 | Ø19 mm reamers kit for 1 st MTP arthrodesis | - Ø19 mm Initial convex reamer - Ø19 mm Initial concave reamer - Pin Ø1.6 L100 mm x 2 |
| KIT-MI22 | Ø22 mm reamers kit for 1 st MTP arthrodesis | - Ø22 mm Initial convex reamer - Ø22 mm Initial concave reamer - Pin Ø1.6 L100 mm x 2 |

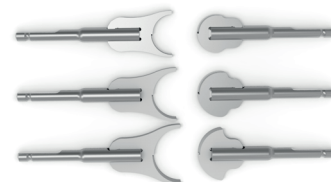

Templates

Sterile templates

| INITIAL F™ - MTP TEMPLATES | |
|----------------------------|--|
| Ref. | Description |
| ANC808 | 1 st MTP arthrodesis plates template |
| ANC846 | 1 st MTP arthrodesis reamers template |

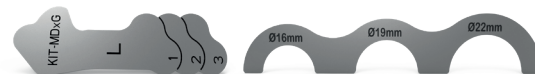


The information presented in this brochure is intended to demonstrate a Newclip Technics product. Always refer to the package insert, product label and/or user instructions before using any Newclip Technics product. Surgeons must always rely on their own clinical judgment when deciding which products and techniques to use with their patients. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Newclip Technics representative if you have questions about the availability of Newclip Technics products in your area.


Ready
when you are!

KIT-MD1D

Example of kit content.

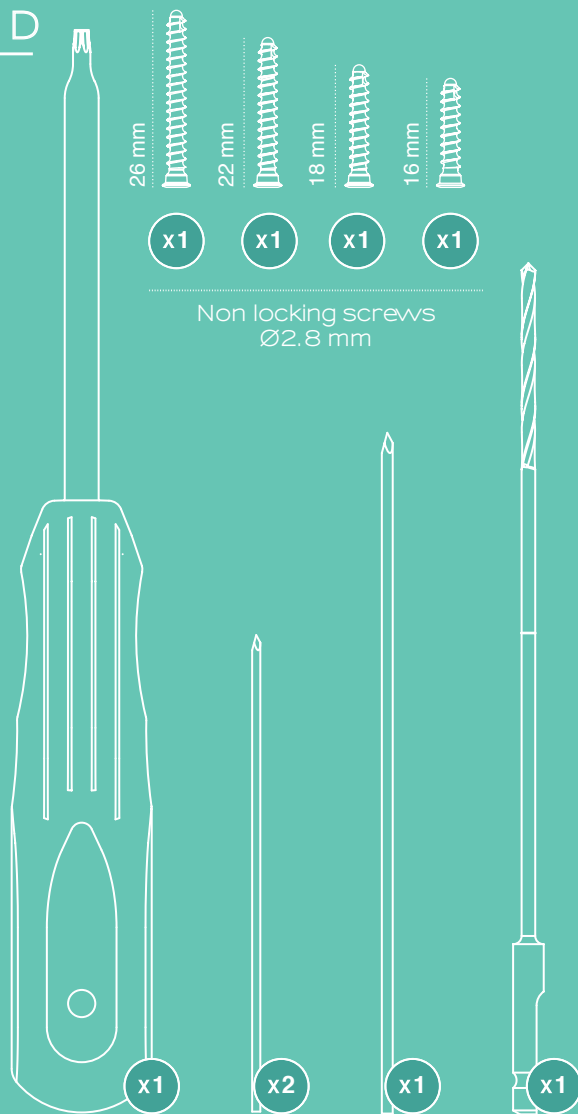

Non locking screws
Ø2.8 mm

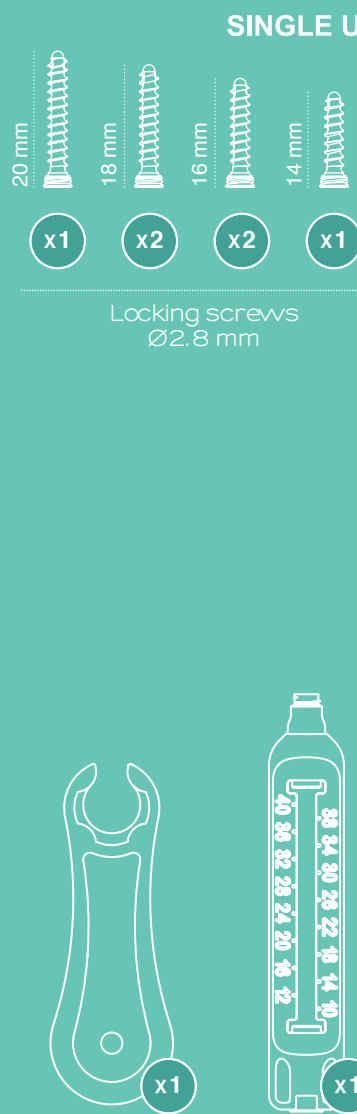

Locking screws
Ø2.8 mm

SINGLE USE KIT

STERILE R




CE
0120

R_x only



Non contractual pictures.

Right foot
1st MTP Arthrodesis
Size 1

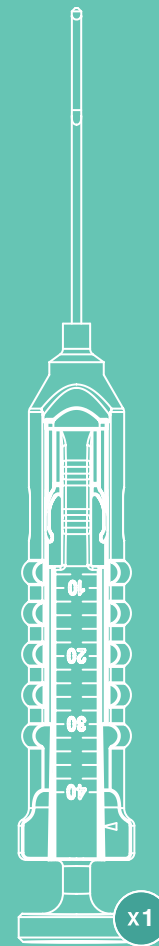

FTM
INITIAL MTP



NEWCLIP TECHNICS
PA de la Lande Saint Martin - 45 rue des Garottières
44115 Haute-Goulaine - France
Tel : +33 (0)2 28 21 37 12 - Fax : +33 (0)2 40 63 68 37
orders@newcliptechnics.com - www.newcliptechnics.com

NEWCLIP USA
642 Larkfield Center
Santa Rosa CA 95403, USA
Phone : + 1 707 230 5078
customerservice@newclipusa.com - www.newclipusa.com

Brochure EN - Initial-F - MTP - Ed.1 - 11/2017 - Medical devices; class IIb - CE0120 SGS UK - US Class: II - Read labeling and instructions before use.