

ACTIV ANKLE MEDIAL MALLEOLUS PLATE



NEW CLIP-TECHNICS

► **Indication :** The implants of the Activ Ankle range are intended for the fixation of fractures, osteotomies and pseudarthroses of the distal and the diaphyseal fibula, for the fixation of fractures of the medial malleolus and for the syndesmotic repair in adults.

Contra-indications:

- Serious vascular deterioration, bone devitalization,
- Pregnancy.
- Acute or chronic local or systemic infections.
- Lack of musculo-cutaneous cover, severe vascular deficiency affecting the concerned area.
- Insufficient bone quality preventing a good fixation of the implants into the bone,
- Muscular deficit, neurological deficiency or behavioral disorders, which could submit the implant to abnormal mechanical strains.
- Allergy to one of the materials used or sensitivity to foreign bodies.
- Serious problems of non-compliance, mental or neurological disorders, failure to follow post-operative care recommendations.
- Unstable physical and/or mental condition.

TECHNICAL FEATURES

MEDIAL MALLEOLUS STANDARD AND HOOK PLATES

- Anatomical implants for the fixation of medial malleolus.
- 2 plates :
 - **Standard** plate with 2 ascendant screws for rigid fixation.
 - **Hook plate** with 1 ascendant screw for a low profile.
- Implant adapted to vertical and horizontal fractures.




STANDARD PLATE

HOOK PLATE

FIXATION

STANDARD PLATE


- Fixation for horizontal fracture

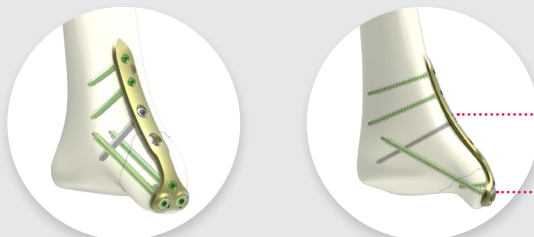

Ramp oblong hole for compression

Ascendant screw through the fracture line

- Fixation for vertical fracture

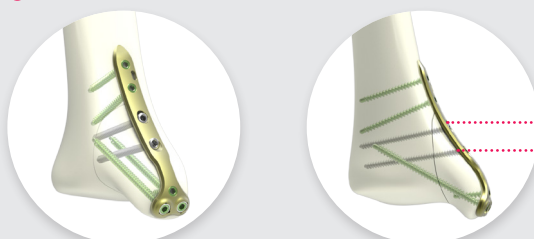


Ramp oblong hole (without compression)


Pre-angulated compression hole

HOOK PLATE

- Fixation for horizontal fracture

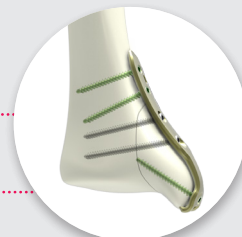


- Fixation for vertical fracture

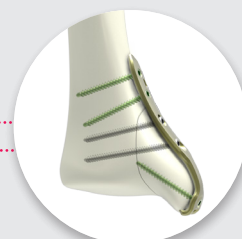

IMPLANTS REFERENCES

PLATES


Ref.	Description
RITSM1	Standard plate for medial malleolus - Size 1
RITSMH1	Hook plate for medial malleolus - Size 1

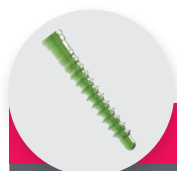

LOCKING SCREWS*

Ref.	Description
SDT2.8L30-STI	Locking screw - Ø2.8 mm - L 30 mm - STERILE
SDT2.8L32-STI	Locking screw - Ø2.8 mm - L 32 mm - STERILE
SDT2.8L34-STI	Locking screw - Ø2.8 mm - L 34 mm - STERILE
SDT2.8L36-STI	Locking screw - Ø2.8 mm - L 36 mm - STERILE
SDT2.8L38-STI	Locking screw - Ø2.8 mm - L 38 mm - STERILE
SDT2.8L40-STI	Locking screw - Ø2.8 mm - L 40 mm - STERILE
SDT2.8L42-STI	Locking screw - Ø2.8 mm - L 42 mm - STERILE
SDT2.8L44-STI	Locking screw - Ø2.8 mm - L 44 mm - STERILE
SDT2.8L46-STI	Locking screw - Ø2.8 mm - L 46 mm - STERILE
SDT2.8L48-STI	Locking screw - Ø2.8 mm - L 48 mm - STERILE
SDT2.8L50-STI	Locking screw - Ø2.8 mm - L 50 mm - STERILE
SDT2.8L55-STI	Locking screw - Ø2.8 mm - L 55 mm - STERILE
SDT2.8L60-STI	Locking screw - Ø2.8 mm - L 60 mm - STERILE

* Green anodized.

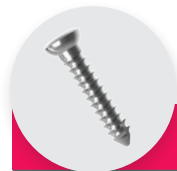

CORTICAL SCREWS*

Ref.	Description
CT2.8L30-ST	Standard cortical screw - Ø2.8 mm - L 30 mm - STERILE
CT2.8L32-ST	Standard cortical screw - Ø2.8 mm - L32 mm - STERILE
CT2.8L34-ST	Standard cortical screw - Ø2.8 mm - L34 mm - STERILE
CT2.8L36-ST	Standard cortical screw - Ø2.8 mm - L36 mm - STERILE
CT2.8L38-ST	Standard cortical screw - Ø2.8 mm - L38 mm - STERILE
CT2.8L40-ST	Standard cortical screw - Ø2.8 mm - L40 mm - STERILE
CT2.8L42-ST	Standard cortical screw - Ø2.8 mm - L42 mm - STERILE
CT2.8L44-ST	Standard cortical screw - Ø2.8 mm - L44 mm - STERILE
CT2.8L46-ST	Standard cortical screw - Ø2.8 mm - L46 mm - STERILE
CT2.8L48-ST	Standard cortical screw - Ø2.8 mm - L48 mm - STERILE
CT2.8L50-ST	Standard cortical screw - Ø2.8 mm - L50 mm - STERILE
CT2.8L55-ST	Standard cortical screw - Ø2.8 mm - L55 mm - STERILE

* Not anodized.

INSTRUMENTATION REFERENCES

INSTRUMENTS

Ref.	Description	Qty
ANC689	Length gauge for Ø2.8 mm screw - Measures 10 - 60 mm	1
ANC690	Ø2.0 mm quick coupling drill bit - L 180 mm	1
ANC691	Ø2.0 mm threaded long guide gauge for Ø2.8 mm screws	1
ANC692	Ø2.0 mm non threaded bent guide gauge for Ø2.8 mm screws	1
ANC268C	Ø2.0 mm threaded guide gauge for Ø2.8 mm screws	1
ANC088C	Ø2.0 mm quick coupling drill bit - L 125 mm	1



NEWCLIP TECHNICS
 PA de la Lande Saint Martin - 45 rue des Garottières
 44115 Haute Goulaine (France)
 Tél. : +33 (0)2 28 21 37 12 - Fax : +33 (0)2 40 63 68 37
 orders@newcliptechnics.com - www.newcliptechnics.com

The information presented in this brochure is intended to demonstrate a NEWCLIP TECHNICS product. Always refer to the package insert, product label and/or user instructions before using any NEWCLIP TECHNICS product. Surgeons must always rely on their own clinical judgment when deciding which products and techniques to use with their patients. Products may not be available in all markets. Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your NEWCLIP TECHNICS representative if you have questions about the availability of NEWCLIP TECHNICS products in your area.